IN240701443V01_GL

713-119V00

EN_IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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③ x 2

⑥ x 2

EN Insert Foot Tube⑥ into Accessories③, making sure the marbles are in the same hole

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EN Insert Accessories② into Accessories③, ensure that the marbles on both sides are in the same hole

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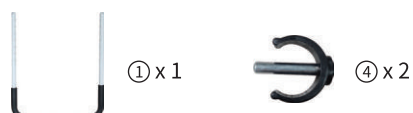
EN Insert Accessories① into Accessories③, make sure the holes on both sides are the same, and then insert C type buckle④

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EN

WARNING

- Please check whether the fittings are complete before installation
- Before use, make sure that the fasteners are properly installed
- When using, be sure to use in a smooth situation
- Please pay attention to safety when using, if necessary, please arrange nursing staff to assist use
- For a period of time each use, please carefully check whether the fastener is loose, the wear of the pad, if necessary, please replace in time
- People with limited mobility or older people should be accompanied by others to prevent secondary injury caused by improper use

MAINTENANCE

- The exterior of should be cleaned periodically with a soft damp cloth. Do not use abrasive cleaners or caustic chemicals.
- Check the parts for any looseness, wear, or damage periodically. If any problems are found, correct them before using them.
- If haven' t used it for a long time, store it in a cool,dry, and safe place.

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If you have any questions, please contact our customer care center, Our contact details are below:

US

 001-877-644-9366

 customerservice@aosom.com

Imported by Aosom LLC
27150 SW Kinsman Rd Wilsonville, OR 97070 USA
MADE IN CHINA

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If you have any questions, please contact our customer care center, Our contact details are below:

CA

 416-792-6088

 customerservice@aosom.ca

Imported by Aosom Canada Inc.
7270 Woodbine Avenue, Unit 307, Markham,
Ontario Canada L3R 4B9
MADE IN CHINA

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If you have any questions, please contact our customer care center, Our contact details are below:

UK

 0044-800-240-4004

 enquiries@mhstar.co.uk

IMPORTER ADDRESS:
MH STAR UK LTD
Unit 27, Perivale Park,
Horsenden lane South
Perivale, UB6 7RH
MADE IN CHINA

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